

Key benefits

- Results in 5-15 minutes in any environment
- Costs 10 times less than drug testing lab
- Minimal training
- Easy and rapid update of drug database
- Hand held and portable, also works from car battery
- Widening list of detected narcotics and chemicals.

Features

- 99.6% accuracy
- Detection of six narcotics
- Two precursors for methamphetamine production
- Differentiation of Cocaine Hydrochloride and Cocaine Base in < 5 min
- Detection limit for controlled substances in street samples = 5
- Fully printable reports.

Technology Overview

- Technology is based on Spectral Fluorescence Signature (SFS) method
- The molecules of a drug of abuse emit a specific light named fluorescence under the excitation of these molecules
- SFS is a fluorescent fingerprint of substances recorded at multi-wavelength excitation of sample
- Sensitivity and selectivity of SFS method considerably exceed the performances of conventional fluorescence technique
- Artificial Neural Networks as expert system recognise the fingerprint of drugs independently of adulterants and diluents in street samples
- Patented technology, unique and revolutionary.

FIND OUT MORE ABOUT NARTEST'S NTX2000 BREAKTHROUGH DRUG PROFILING SYSTEM

If you're ready to take charge with fast, accurate testing of seized drug samples, you're ready to try NarTest NTX2000.

To find out more about how your organisation can put the NarTest NTX 2000 Drug Analyser to work, contact us today to arrange a demonstration.



Arrivals hall



Customs


Douane



Fast, accurate, real-time Drugs Identification

NarTest NTX2000

Drug Testing

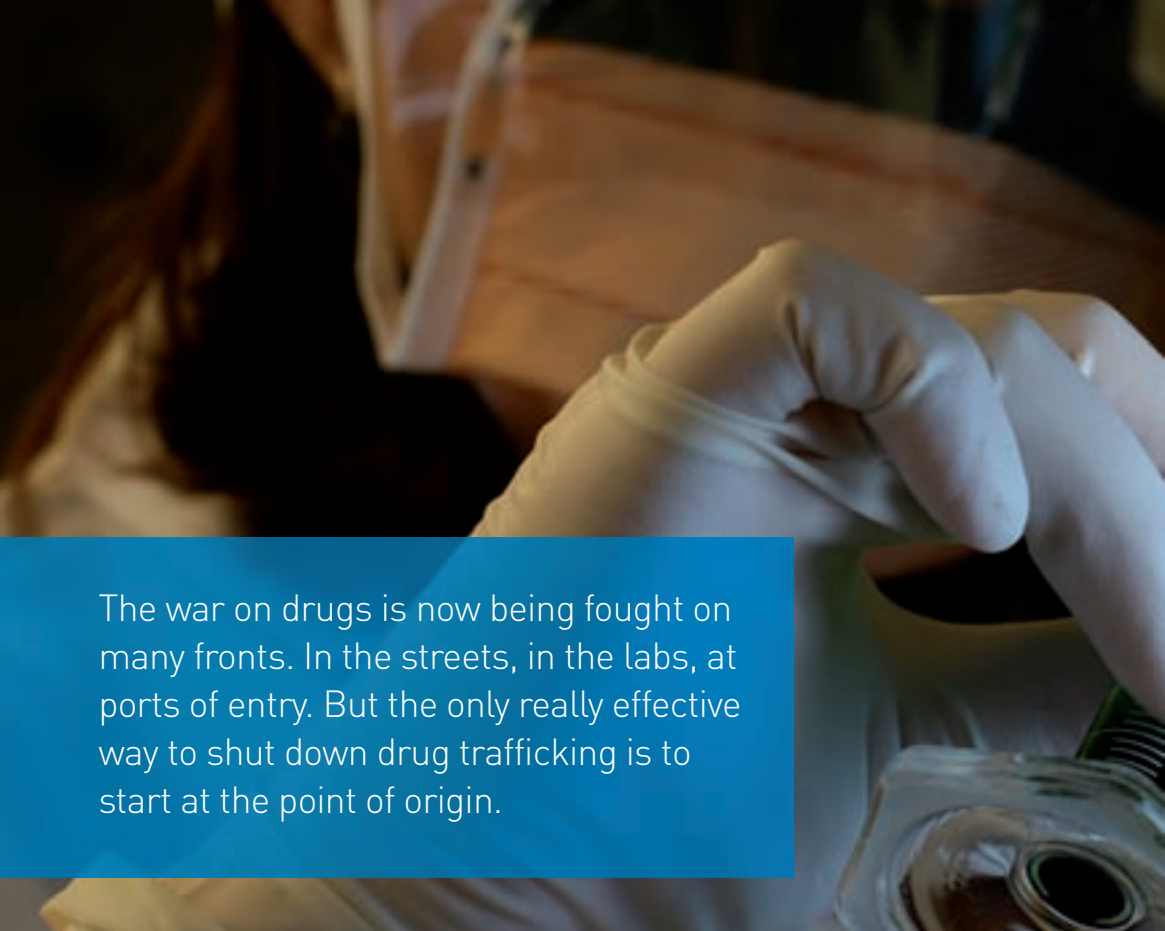

The war on drugs is now being fought on many fronts. In the streets, in the labs, at ports of entry. But the only really effective way to shut down drug trafficking is to start at the point of origin.

NarTest NTX2000 can shorten the time between arrest and the court room

Drug testing is a difficult process for all concerned, primarily because few laboratory facilities are located nearby. In addition, almost every laboratory faces a large number of samples awaiting testing. Inevitably, hard-core drug dealers can, and will, work their trade in the interim.

Even smaller dealers can easily rack up hundreds of thousands of pounds in illicit drug trade while they're still active. Most pushers know how to take advantage of any waiting period for testing, well aware that a high-speed, reliable drug-testing system has not been available to law enforcement.

Hard evidence. Low cost.

NarTest NTX2000 uses a completely self-contained instrument. Its compact size and user-friendly system mean it's easily deployed in police stations and local or mobile lab facilities. Full training and technical support is available. So that local personnel can confidently operate the system.

With its demonstrated speed and unblemished record for accuracy, NarTest NTX2000 means fast decisions about drugs.

It's one thing to seize illegal drugs.
Now you can find out where they came from.

Proof positive, every time.



Fast, effective, proven

Thanks to proven Spectral Fluorescence Signature (SFS) technology, the NarTest NTX2000 system can analyse evidence in a matter of minutes. Cocaine, heroin, methamphetamine, marijuana and ecstasy - accounting for 84% of the most frequently-identified illegal drugs - can be accurately profiled in the field and lead authorities along the entire route from manufacture to seizure.

In approximately ten minutes, samples prepared by easily-trained, non-scientific local personnel can be analysed using this patented technology. Using the NarTest NTX2000 substance database, even complex mixtures at low drug concentration can be readily matched, and identified in STGT real time. Not only does the instrument alert the operator immediately to the presence of an illegal drug, the technology then generates a date/time-stamped identification report for the case file.